Quick Start Guide 1

The program buttons and home button have 2 button settings; 1 Momentary (press & release) or 2 Maintain (press & hold).
For instructions on changing the button settings consult the user manual.


		To Save a Position	To Move to a Position
1 Momentary Setting Factory Default	PROGRAM BUTTONS P1 P2 P3 P4	Hold any P1 2 3 4 5 sec	Press & Release any P1 2 3 4
	PROGRAMMABLE HOME	Hold 10 sec	Press & Release
	ADA HOME	Factory Setting (not changeable)	Double Press → Release
2 Maintain Setting	PROGRAM BUTTONS P1 P2 P3 P4	→ → Hold any P1 2 3 4 5 sec	Press & Hold any P1 2 3 4
	PROGRAMMABLE HOME	→ → Hold 10 sec	Press & Hold
	ADA HOME	Factory Setting (not changeable)	Press → Release → Press & Hold

MTI, Inc. • 3655 West Ninigret Drive, Salt Lake City, Utah 84104 USA • (877) 875-4990 • +1 (801) 875-4990 • technical.support@mti.net • www.mti.net


To prevent unintentional operation of the chair use the following commands:

SAFETY LOCKOUT	To Lock	To Unlock
		

The swivel button operates differently depending on the button setting as shown. See User Manual to change the button setting.

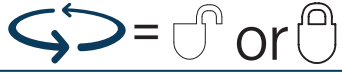SWIVEL	Auto-Lock Setting	Lock/Unlock Setting
		

INDICATOR LIGHTS

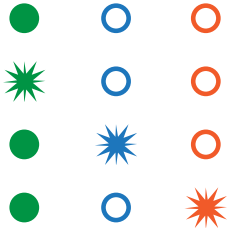


Chair is operating normally

Chair is hibernating

Chair is moving

Swivel is unlocked

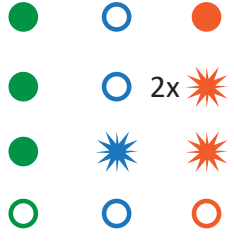

Safety Lockout is enabled

2x

With accompanying beep, chair has reached boundary of safe zone or excessive load

Chair is moving to ADA Compliant Home

Chair is off

WARNING/AVERTISSEMENT

Read the user manual before installation and use.
Lisez le manuel d'utilisation avant l'installation et l'utilisation.

WARNING/AVERTISSEMENT

Ensure equipment, patient, and other objects are in safe condition before operating the equipment.
Assurer l'appareil, patient, et d'autre objets sont dans un état sûr avant d'utiliser l'appareil.

© MTI All rights reserved. MTI reserves the right to make any product changes without notice. The MTI brand and Bear logo design are U.S. registered trademarks of MTI, Inc.