Quick Start Guide 2 - Button Operation

The program buttons and home button have 2 button settings; 1 Momentary (press & release) or 2 Maintain (press & hold).
For instructions on changing the button settings consult the user manual.

1 Momentary Button Setting (Press & Release) - Factory Default		
	To Save a Position	To Move to a Position
PROGRAM BUTTONS P1 P2 P3 P4	Hold any P1 2 3 4 5 sec	Press & Release any P1 2 3 4
PROGRAMMABLE HOME	Hold 10 sec	Press & Release
ADA HOME	Factory Setting (not changeable)	Double Press → Release

2 Maintain Button Setting (Press & Hold)		
	To Save a Position	To Move to a Position
PROGRAM BUTTONS P1 P2 P3 P4	→ → Hold any P1 2 3 4 5 sec	Press & Hold any P1 2 3 4
PROGRAMMABLE HOME	→ → Hold 10 sec	Press & Hold
ADA HOME	Factory Setting (not changeable)	Press → Release → Press & Hold


To prevent unintentional operation of the chair use the following commands:		
SAFETY LOCKOUT	To Lock → → →	To Unlock → → →

INDICATOR LIGHTS

- Chair is operating normally
- Chair is hibernating
- Chair is moving
- Swivel button was pressed
- Safety Lockout is enabled
- With accompanying beep, chair has reached boundary of safe zone or excessive load (2x)
- Chair is moving to ADA Compliant Home
- Chair is off

BATTERY OPERATION

The battery must be firmly inserted with the label pull tab facing out. If the chair is not plugged in to a wall outlet, the bed will play a start-up tune when properly inserted.

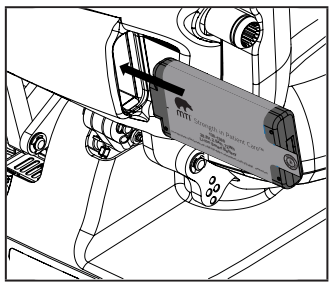

WARNING/AVERTISSEMENT

Read the user manual before installation and use.
Lisez le manuel d'utilisation avant l'installation et l'utilisation.

WARNING/AVERTISSEMENT


Ensure equipment, patient, and other objects are in safe condition before operating the equipment.
Assurer l'appareil, patient, et d'autre objets sont dans un état sûr avant d'utiliser l'appareil.

Quick Start Guide 2 - Steering Options

The 840 Stretcher Bed has 3 different steering options: 1 Rear Steer (Model 840-003-22), 2 Front Steer (Model 840-003-23), and 3 Selective Steer (Model 840-003-24). Check the last two digits of the model number on the serial label of the stretcher bed to identify the steering option. Each steering option has 3-4 modes. Use the guide below to set the steering bars in the desired mode.

1 Rear Steer


The Rear Steer Option utilizes the front caster on the patient's left side as the Steering Caster and the other casters are non-steering casters.

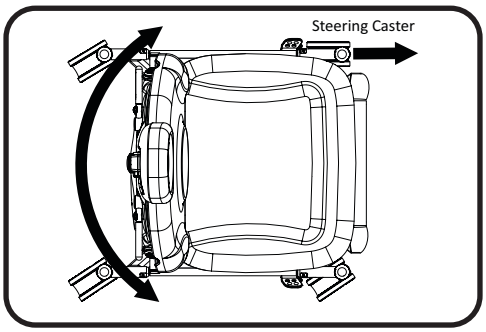

Mode	Left Side	Right Side
Steering Mode Allows rear of bed to swivel		
Free Mode Allows all casters to swivel		
Lock Mode Locks all casters		

2 Front Steer


The Front Steer Option utilizes the rear caster on the patient's right side as the Steering Caster and the other casters are non-steering casters.

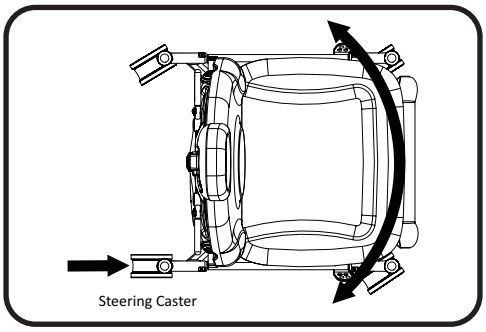

Mode	Left Side	Right Side
Steering Mode Allows front of bed to swivel		
Free Mode Allows all casters to swivel		
Lock Mode Locks all casters		

3 Selective Steer

The Selective Steer Option utilizes the rear caster on the patient's right side and the front caster on the patient's left side as Steering Casters and the other casters are non-steering casters.

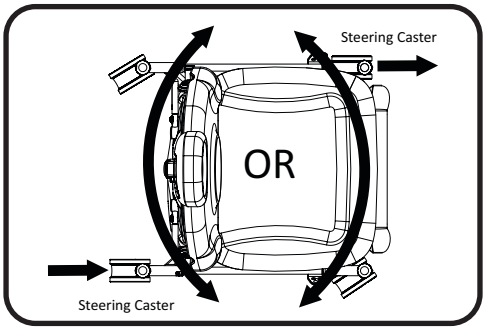

Mode	Left Side	Right Side
Rear Steering Mode Allows rear of bed to swivel		
Front Steering Mode Allows front of bed to swivel		
Free Mode Allows all casters to swivel		
Lock Mode Locks all casters		

⚠ *WARNING/AVERTISSEMENT*

If the stretcher bed has the Selective Steer Option, do not put both the left-hand and right-hand locking bar in the STEER Position. This will prevent the stretcher bed from turning.

Si la civière mobile a l'option de Conduire Sélective, ne mettez pas la barre de verrouillage à droite et à gauche dans la position de STEER (CONDUIRE) à la fois. Cela interdire la civière mobile de tourner.

MTI, Inc. • 3655 West Nanigret Drive, Salt Lake City, Utah 84104 USA
(877) 875-4990 • +1 (801) 875-4990
technical.support@mti.net • www.mti.net

© MTI All rights reserved. MTI reserves the right to make any product changes without notice. The MTI brand and Bear logo design are U.S. registered trademarks of MTI, Inc.

130-3274 Rev B